

ultrapure water tailored for trace elemental analysis

Milli-Q[®] IQ Element

Water Purification and Dispensing Unit



Milli-Q



The Life Science Business of Merck operates as MilliporeSigma in the US and Canada.



ultrapure water without compromise

Milli-Q® IQ Element Water Purification and Dispensing Unit

Ultrapure water suitable for the most stringent trace elemental analyses

Don't let trace contaminants interfere with your sensitive analyses.

The Milli-Q[®] IQ Element unit, combined with a Milli-Q[®] IQ 7 series water purification system, delivers analytical-grade ultrapure water that is suitable for trace and ultra-trace elemental analyses, including **ICP-MS**, **GF-AAS** and **trace IC**.

Using fresh ultrapure water from a Milli-Q[®] IQ 7000 or Milli-Q[®] IQ 7003/05/10/15 system, the Milli-Q[®] IQ Element unit purifies even further. The water dispensed at your point of use is confirmed to contain **extremely low levels of elemental contaminants, from single ppt to sub-ppt detection levels**^{*}. Independent laboratories specialized in ultra-trace elemental analyses have verified the unit's water quality.

Designed to fit into your trace analysis workflow

Easy to Integrate

The compact unit is designed for seamless, contaminant-free installation in your cleanroom environment or laminar flow hood.

Easy to Use

A touchscreen lets you continuously view essential quality parameters and, in a few clicks, you can print a dispense report or program your desired dispense volume.

Easy to Avoid Contamination

There's no need to touch the unit while working; a footswitch allows for hands-free dispensing at your point of use.

Easy to Maintain

All purification cartridges have been designed to be effortlessly replaced. The procedure can be performed without the intervention of a field service engineer.

Easy Data Management

Never lose track of your water quality. An intuitive data management system lets you monitor, store and rapidly retrieve water quality data in a few clicks—from a single dispense to a complete history.



Designed to produce and maintain the purity of high-quality ultrapure water

Removes trace ionic contaminants from feed ultrapure water

Milli-Q® IQ 7 series water purification system Delivers consistently high-quality ultrapure water




IPAK Quanta® ICP polishing cartridge Removes trace ions



0.1 μm Optimizer LW™ final filter Removes trace particulates Product ultrapure water is suitable for sensitive ICP-MS applications.



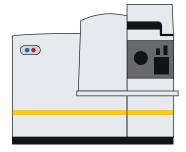

Safeguards against the introduction of contaminants

The Milli-Q[®] IQ Element unit not only further purifies ultrapure water down to trace (ppt) and ultra-trace (sub-ppt) levels, but its design protects the water from the introduction of contaminants from the environment.

- All components used for water production are made from **selected low-extractable materials**
- **Footswitch and dispenser** provide hands-free water delivery to reduce risk of contamination from your surroundings while you are working
- Touchscreen display allows for water quality monitoring at a glance







Intuitive touchscreen allows for easy control, monitoring and maintenance






IQ Element

Quality Monitoring



Volumetric Dispensing

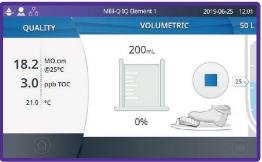

Data Management



Maintenance Wizards



Easy integration into your laboratory space

The compact Milli-Q[®] IQ Element unit can be easily installed inline of a Milli-Q[®] IQ 7 series water purification system. Its simple dispensing unit can be placed directly at your point of use, in a clean and controlled environment, with no risk of contamination.



technical appendix

Ultrapure (Type 1) Water Specifications

Resistivity @ 25°C 18.	3.2 MΩ.cm
Total Organic Carbon (TOC) \leq 5	5 ppb
Flow rate up	to 1.5 L/min

1. Under standard operating conditions. For more information, please refer to the user manual of the Milli-Q® IQ 7000 or Milli-Q® IQ 7003/05/10/15 water purification systems.

Dimensions & weights

Parameter	Value			
Dimensions (H x W x D)	67.1 x 21.1 x 27.0 cm 26.4 x 8.3 x 10.6 in.			
Dry weight	7.5 kg (16.5 lb)			
Shipping weight	10.3 kg (22.7 lb)			
Operating weight	9.1 kg (20.1 lb)			

Ordering information

Description	Catalog number
Milli-Q [®] IQ Element unit	ZIQELEMT0
Milli-Q [®] IQ Element consumables kit	IPAKICPK1

ICP-MS analysis of ultrapure water from Milli-Q® IQ Element unit

Excerpt from ICP-MS analysis with high purity water obtained from a Milli-Q[®] IQ Element purification unit connected to a Milli-Q[®] IQ 7005 water purification system. Additional results as well as detailed experimental methods are available in the Milli-Q[®] IQ Element Data Sheet.

Isotope	Element	Sample (ng/L)	DL (ng/L)	Isotope	Element	Sample (ng/L)	DL (ng/L)
7	Lithium (Li) ¹	< DL	0.04	115	Indium (In) ²	< DL	0.01
9	Beryllium (Be) ²	< DL	0.20	118	Tin (Sn) ²	< DL	0.15
11	Boron (B) ²	< DL	0.50	121	Antimony (Sb) ¹	< DL	0.02
23	Sodium (Na)1	0.68	0.11	126	Tellerium (Te) ²	0.08	0.07
24	Magnesium (Mg) ¹	0.01	0.01	133	Caesium (Cs) ¹	0.01	0.00
27	Aluminium (Al) ¹	0.07	0.04	138	Barium (Ba) ¹	< DL	0.05
28	Silicon (Si) ¹	198.65*	4.98	139	Lanthanium (La) ²	< DL	0.02
39	Potassium (K) ¹	0.54	0.16	140	Cerium (Ce) ²	< DL	0.03
40	Calcium (Ca) ²	< DL	0.29	141	Praseodymium (Pr) ²	< DL	0.02
45	Scandium (Sc) ²	0.59	0.53	146	Neodenyum (Nd) ²	< DL	0.08
47	Titanium (Ti) ¹	0.61	0.51	147	Samarium (Sm) ²	< DL	0.13
51	Vanadium (V) ¹	0.03	0.01	153	Europium (Eu) ²	< DL	0.04
52	Chromium (Cr) ¹	0.08	0.02	157	Gadolinium (Gd) ²	< DL	0.13
55	Manganese (Mn) ¹	0.01	0.02	159	Terbium (Tb) ²	< DL	0.02
56	Iron (Fe) ¹	< DL	0.50	163	Dysprosium (Dy) ²	< DL	0.07
59	Cobalt (Co) ¹	< DL	0.01	165	Holmium (Ho) ²	< DL	0.02
60	Nickel (Ni) ¹	< DL	0.16	166	Erbium (Er) ²	< DL	0.11
63	Copper (Cu) ¹	< DL	0.04	169	Thulium (Tm) ²	< DL	0.03
66	Zinc (Zn) ¹	< DL	0.48	172	Ytterbium (Yb) ²	< DL	0.09
70	Germanium (Ge) ²	< DL	0.10	175	Lutetium (Lu) ²	< DL	0.02
71	Gallium (Ga) ²	< DL	0.13	178	Hafnium (Hf) ²	< DL	0.11
75	Arsenic (As) ¹	0.06	0.04	181	Tantalum (Ta) ²	< DL	0.03
78	Selenium (Se) ²	< DL	0.57	182	Tungsten (W) ¹	< DL	0.07
85	Rubidium (Rb) ¹	< DL	0.03	185	Rhenium (Re) ²	< DL	0.09
88	Strontium (Sr) ²	< DL	0.02	189	Osmium (Os) ²	< DL	0.14
89	Yttrium (Y) ²	< DL	0.02	193	Iridium (Ir) ²	< DL	0.05
90	Zirconium (Zr) ²	< DL	0.05	195	Platinum (Pt) ²	0.18	0.16
93	Niobium (Nb) ²	< DL	0.03	197	Gold (Au) ²	< DL	0.43
95	Molybdenium (Mo) ¹	< DL	0.10	202	Mercury (Hg) ²	5.1	1.52
101	Ruthenium (Ru) ²	0.42	0.20	205	Thallium (Tl) ²	<dl< td=""><td>0.05</td></dl<>	0.05
103	Rhodium (Rh) ²	< DL	0.01	208	Lead (Pb) ¹	< DL	0.08
105	Palladium (Pd) ²	< DL	0.34	209	Bismuth (Bi) ²	< DL	0.06
107	Silver (Ag) ²	0.40	0.15	232	Thorium (Th) ²	<dl< td=""><td>0.04</td></dl<>	0.04
111	Cadmium (Cd) ¹	< DL	0.08	238	Uranium (U) ²	<dl< td=""><td>0.04</td></dl<>	0.04

DL: Detection limit

1. Data obtained courtesy of Agilent Technologies, Tokyo, Japan. © Agilent Technologies, Inc. Reproduced with Permission, Courtesy of Agilent Technologies, Inc.

2. Data obtained courtesy of UT2A, Pau, France.

* Si is known to be difficult to measure by ICP-MS. When measured by GF-AAS, concentration was < DL (0.5 ppb).

Milli-Q ® Lab Water Solutions

 $\textcircled{\sc c}$ 2019 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved. Merck, the vibrant M, Milli-Q, IPAK Quanta and Optimizer LW are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

MK_BR4224EN

